



Auto Lensmeter LM-600

Serving Your Vision



Auto Lensmeter LM-600

Wavefront technology with Hartmann sensor (117 points)

Using 117 point Hartmann sensor wavefront technology in the LM-600, measurement accuracy of vertically narrow progressive lenses has markedly improved.

Simultaneous measurement of UV/blue light and power


By emitting measurement light, UV light (375 nm) and blue light (465 nm) simultaneously, it is possible to check UV/blue light transmittance as well as conduct power measurements.

Lens mark recognition support

Displaying the grid patterns on the screen will help recognition of lens marks. Grid patterns are selectable.


LAN and RS-232C connection

Output the measurement data from each port of LAN and RS-232C.

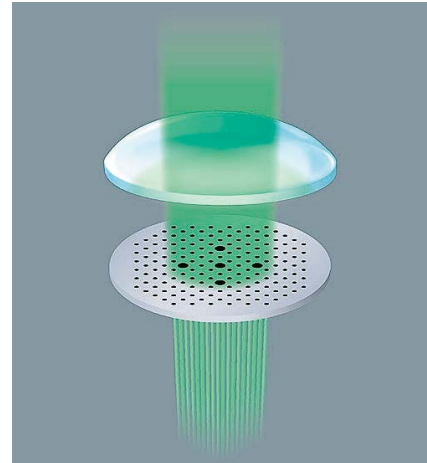

Specifications

Measurement range		Measurement object	Spectacle lens, Contact lens, Optical lens
Spherical power(SPH)	±25 D	Transmittance of UV	The peak of the wavelength is 375 nm
Cylindrical power(CYL)	±10 D	Transmittance of blue light	The peak of the wavelength is 465 nm
Axial angle(Axis)	0-180°	Display	7-inch-wide color TFT with touch panel
Additional power(ADD)	-2 to +10 D	Printer	Thermal printer, Paper width 58 mm
Prism power	0-15 Δ	External communication	RS-232C, USB3.0, Ethernet
Increment		Dimensions and power source	
Diopter	0.01, 0.06, 0.12, 0.25 D	Dimensions (when LCD is tilted)	188 (W) × 240 (D) × 430 (H) mm
Prism	0.01, 0.06, 0.12, 0.25 Δ	Weight	Approx. 5.5 kg
Wavelength	535 nm	Voltage	AC 100-240 V
Diameter of the lens	Ø20-120 mm (more than Ø5 mm for CL)	Frequency	50/60 Hz
		Power consumption	40-50 VA

 To ensure safe and correct handling please read the user manual before using.

- Description and appearance as detailed in this brochure may be subject to change as improvements are made to products.
- Colours as they appear in brochure photographs may be slightly different to actual products due to lighting conditions when photographs were taken, or print colours.

 For The Americas, Asia-Pacific & Middle East
TAKAGI SEIKO CO., LTD.
330-2 Iwafune, Nakano-shi, Nagano-ken, 383-8585, Japan
TEL : +81(0)269-22-4511(Switchboard) URL : <https://www.takagi-j.com>



For Europe & Africa
Takagi Ophthalmic Instruments Europe Ltd
Citylabs 1.0, Nelson Street, Manchester, M13 9NQ, UK
TEL : +44 (0)161 273 6330 URL : <https://www.takagieurope.com>

**CEUK
CA**

Your local distributor